Compact laboratory balance KERN PFB











Quick-displaying precision balance with user-friendly operation modes and EC type approval

Features

- Easy to use: All primary functions have their own key on the keypad
- Compact size, practical for small spaces
- Capacity display: A bar lights up to show how much of the weighing range is still available
- Level indicator and levelling feet for precise levelling of the scale, fitted as standard, to give the most accurate weighing results
- Draught shield standard for models with weighing plate size A and B. Weighing space WxDxH 158x143x64 mm

• Separate second display standard for models PFB 600-1M and PFB 6000-0M. Suitable for direct sales to the public. The verification value is the same as the readout. Large backlit LCD display, digit height 16 mm. Interface cable to the balance, cable length approx. 1 m, standard.

Dimensions WxDxH 250x180x105 mm

Technical data

- Backlit LCD display, digit height 15 mm
- Weighing plate dimensions, stainless steel, A Ø 80 mm,
- **B** Ø 120 mm,
- WxD 155x145 mm, see enlarged picture

- Overall dimensions (without draught shield) WxDxH 200x190x87 mm
- Net weight approx. 2 kg
- \bullet Permissible ambient temperature 5 °C / 35 °C

Accessories

- Separate second display, ideal for training and demonstration purposes in laboratories or industry, not permitted for direct sales to the public, KERN PFB-A08
- Bluetooth data interface for wireless data transfer to PC, must be ordered at purchase, KERN PFB-A10
- Suitable printers see page 177 ff.

STANDARD































Model	Weighing	Readout	Verification	Repro-	Linearity	Weighing	Options			
	range		value	ducibility		plate	Verification		DAkkS Calibr. Certificate	
	[Max]	[d]	[e]				MO		DKD	
KERN	g	g	g	g	g		KERN		KERN	
PFB 120-3	120	0,001	-	0,001	± 0,003	Α	-	-	963-127	
PFB 200-3	200	0,001	-	0,001	± 0,003	А	-	ı	963-127	
PFB 300-3	300	0,001	-	0,002	± 0,003	А	-	-	963-127	
PFB 1200-2	1200	0,01	-	0,01	± 0,03	В	-	-	963-127	
PFB 2000-2	2000	0,01	-	0,01	± 0,03	В	-	-	963-127	
PFB 3000-2	3000	0,01	-	0,01	± 0,03	В	-	-	963-127	
PFB 6K0.05	6000	0,05	-	0,05	± 0,15	С	-	ı	963-128	
PFB 6000-1	6000	0,1	-	0,1	± 0,3	С	-	-	963-128	

Note: For applications that require verification, please order verification at the same time, initial verification at a later date is not possible. Verification at the factory, we need to know the full address of the location of use.

PFB 600-2M	600	0,01	0,1	0,01	± 0,02	В	965-216	963-127	
PFB 600-1M	600	0,1	0,1	0,1	± 0,1	В	965-216	963-127	
PFB 6000-1M	6000	0,1	1	0,1	± 0,3	С	965-217	963-128	
PFB 6000-0M 🔳	6000	1	1	1	± 1	C	965-217	963-128	

KERN Pictograms:



Internal adjusting: Quick setting up of the balance's accuracy with internal adjusting weight (motordriven).



Adjusting program CAL: For quick setting up of the balance's accuracy. External adjusting weight required.



Memory: Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.



Alibi memory: Electronic archiving of weighing results, complying with the 2009/23/EC standard.



Data interface RS-232: To connect the balance to a printer, PC or network.



RS-485 data interface: To connect the balance to a printer, PC or other peripherals. High tolerance against electromagnetic disturbance.



USB data interface: To connect the balance to a printer, PC or other peripherals.



Bluetooth* data interface: To transfer data from the balance to a printer, PC or other peripherals.



WLAN data interface: To transfer data from the balance to a printer, PC or other peripherals.



Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.



Interface for second balance: For direct connection of a second balance.



Network interface: For connecting the scale to an Ethernet network. With KERN products you can use a universal RS-232/LAN converter.



Wireless data transfer: between the weighing unit and the evaluation unit using an integrated radio module.



GLP/ISO log: The balance displays the weight, date and time, regardless of a printer



GLP/ISO log: With weight, date and time. Only with KERN printers.



Piece counting: Reference quantities selectable. Display can be switched from piece to weight.



Recipe level A: Separate memory for the weight of the tare container and the recipe ingredients (net total).



Recipe level B: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display.



Recipe level C: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display, adjustment of recipe when dosages are exceeded, multiplier function, barcode.



Totalising level A: The weights of similar items can be added together and the total can be printed out.



Totalising level C: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display, adjustment of recipe when dosages are exceeded, multiplier function, barcode recognition.




Percentage determination: Determining the deviation in % from the target value (100 %).



model. Please refer to KFRN's website for more details. Weighing with tolerance range: Upper and lower limiting values can be programmed indivi-

Weighing units: Can be switched to e.g. non-

metric units at the touch of a key. See balance

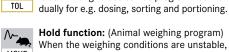

MOVE

Hold function: (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average



Protection against dust and water splashes IPxx: The type of protection is shown in the pictogram.



ATEX explosion protection: Suitable for use in hazardous industrial environments, in which there is explosion danger. The ATEX marking is specified for each device.



Stainless steel: The balance is protected against corrosion.



Suspended weighing: Load support with hook on the underside of the balance.



Battery operation: Ready for battery operation. The battery type is specified for each device.



Rechargeable battery pack:

Rechargeable set.



Universal mains adapter: with universal input and optional input socket adapters for

A) EU, GB B) EU, GB, CH, USA C) EU, GB, CH, USA, AUS



Mains adapter: 230V/50Hz in standard version for EU. On request GB, USA or AUS version available.



Power supply: Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request.



Weighing principle: Strain gauge Electrical resistor on an elastic deforming body.



Weighing principle: Tuning fork A resonating body is electromagnetically excited, causing it to oscillate.



Weighing principle: Electromagnetic force compensation Coil inside a permanent magnet. For the most accurate weighings.



Weighing principle: Single cell technology Advanced version of the force compensation principle with the highest level of precision.



Verification possible:

The time required for verification is specified in the pictogram.



DAkkS calibration possible (DKD): The time required for DAkkS calibration is shown in days in the pictogram.



Package shipment: The time required for internal shipping preparations is shown in days in the pictogram.



Pallet shipment: The time required for internal shipping preparations is shown in days in the pictogram.



Warranty: The warranty period is shown in the pictogram.

KERN – Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg - 2000 kg. In combination with a DAkkS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAkkS calibration laboratory today is one of the most modern and best-equipped DAkkS calibration laboratories for balances, test weights and forcemeasurement in Europe.

Thanks to the high level of automation, we can carry out DAkkS calibration of

balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

Range of services:

- DAkkS calibration of balances with a maximum load of up to 50 t
- DAkkS calibration of weights in the range of 1 mg 2500 kg
- · Database supported management of checking equipment and reminder service Calibration of force-measuring devices
- DAkkS calibration certificates in the following languages D, GB, F, I, E, NL, PL

Your KERN specialist dealer: